

Rabbit anti-human PTEN Polyclonal Antibody

Catalog Number: R16056P

General Information

Immunogen	A synthetic peptide of human PTEN
IgG type	IgG
Clonality	Polyclonal
Specificity	human PTEN
Applications & dilution	WB 1:500 - 1:2000 IHC 1:50 - 1:200
Formulation	PBS with 0.02% sodium azide, 50% glycerol, pH7.3.
Purity	≥95% purity by SDS-PAGE
Storage	Store at -20°C. Avoid freeze / thaw cycles.

Abbreviation:

ELISA: Enzyme-linked immunosorbent assay; ITA: immunoturbidimetric assay; IP: immunoprecipitation; IHC: immunohistochemistry; IF: immunofluorescence. WB: western blot; FC: flowcytometry

Background

PTEN (phosphatase and tensin homologue deleted on chromosome ten), also referred to as MMAC (mutated in multiple advanced cancers) phosphatase, is a tumor suppressor implicated in a wide variety of human cancers. PTEN encodes a 403 amino acid polypeptide originally described as a dual-specificity protein phosphatase. The main substrates of PTEN are inositol phospholipids generated by the activation of the phosphoinositide 3-kinase (PI3K). PTEN is a major negative regulator of the PI3K/Akt signaling pathway. PTEN possesses a carboxy-terminal, noncatalytic regulatory domain with three phosphorylation sites (Ser380, Thr382, and Thr383) that regulate PTEN stability and may affect its biological activity. PTEN regulates p53 protein levels and activity and is involved in G protein-coupled signaling during chemotaxis.


For research use only

Preparation

Polyclonal antibody is produced by immunizing rabbit with a synthetic peptide of human PTEN and purified using protein A resin.

Applications

Western blot

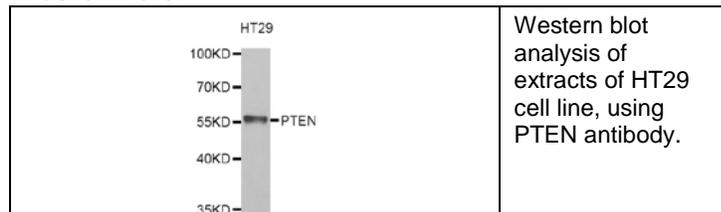

Storage

This antibody is shipped at 4 °C. This product is stable for 12 months from date of receipt when stored at -20 °C to -70 °C. Avoid freeze/thaw cycles.

Hazard/Biohazard

This antibody contains 0.02% sodium azide as preservative. Please handle and dispose the product properly. No known biohazard is associated with this product.